



### NaviTip™ 31 ga with Sideport Irrigator Tips

- Double sideports deliver irrigants safely, minimizing the possibility of chemicals being expressed past the apex
- The world's smallest cannula navigates the most intricate canal spaces

Designed for: ChlorCid, ChlorCid V, ChlorCid Surf, Ultradent EDTA 18% Solution, and Citric Acid 20%.



### NaviTip™ FX™ Tips

- One-of-a-kind brush cleans, scrubs, and irrigates simultaneously
- Rigid 30 ga cannula



Listed as an "EXCELLENT" product by a prominent independent research institute.<sup>1</sup>

Designed for: Ultradent Citric Acid 20% Solution.

<sup>1</sup> Clinical Research Associates Newsletter. Volume 29, Issue 1, January 2005.

31 ga (20pk)  
5121 - 21 mm  
5123 - 27 mm

31 ga (50pk)  
5122 - 21 mm  
5124 - 27 mm

1452 - 17 mm 20pk

1454 - 25 mm 20pk



### NaviTip™ Tips

- Provide controlled delivery to the apex
- Flexible, stainless steel cannula easily navigates curved canals
- Available in 17, 21, 25, and 27 mm lengths and 29 or 30 ga cannulae

29 ga designed for: MTAFlow, EndoREZ, and UltraCal XS.  
30 ga designed for: File-Eze, Ultradent EDTA 18% Solution, and Consepis.



### Endo-Eze™ Tips

- Work for a variety of endodontic procedures
- Come with flexible, strong cannulae in a variety of sizes
- Bend easily

Designed for: luting materials and air/water delivery.  
Use with: TriAway Adapter, PermaFlo DC (20 ga), and other Ultradent syringes.

|               | 29 ga (20pk) | 30 ga (20pk) | 29 ga (50pk) | 30 ga (50pk) |
|---------------|--------------|--------------|--------------|--------------|
| 17 mm, White  | 5112         | 1249         | 1378         | 1421         |
| 21 mm, Yellow | 5113         | 1349         | 1374         | 1422         |
| 25 mm, Blue   | 5114         | 1250         | 1376         | 1423         |
| 27 mm, Green  | 5115         | 1354         | 1377         | 1424         |
| Assorted      | 5116         | 1351         | 1379         | 3319         |


20pk  
345  
346  
347  
348

100pk  
1428  
1429  
1430  
1431

Gauge - Diameter  
18 ga - 0.0490";1.25 mm  
19 ga - 0.0415";1.06 mm  
20 ga - 0.0350";0.90 mm  
22 ga - 0.0275";0.70 mm



Quantity discounts available at 3 and 5 packs.

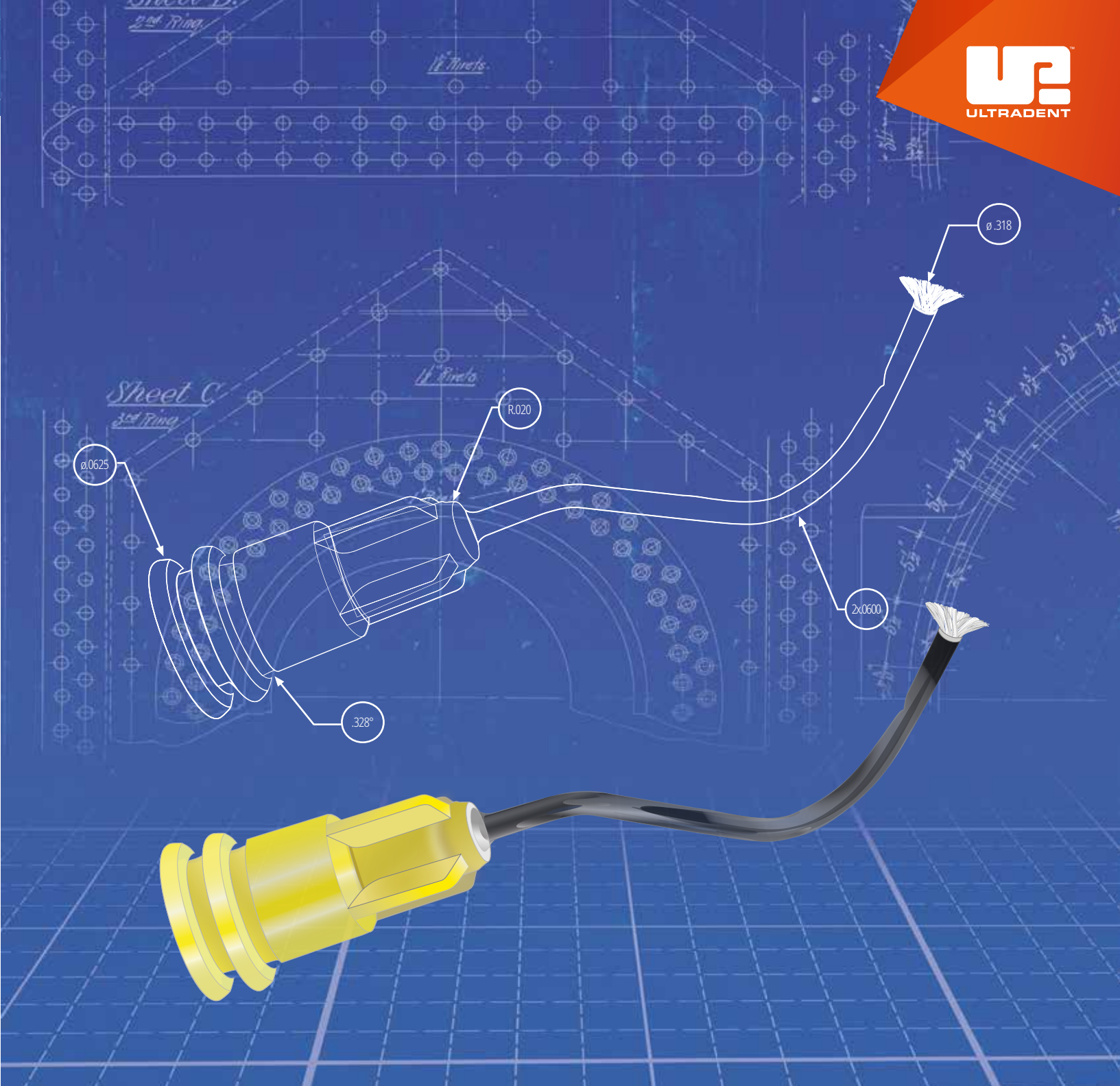

ULTRADENT TIPS

**DESIGNED TO DELIVER**

The chemistries you use are different. Some are chemically activated, needing to be mixed immediately before delivery. Others have varying viscosities. Some work in pits and fissures, some inside canals, and some on enamel surfaces. Each chemistry you use is designed for a specific purpose. Shouldn't the same be true for your tips? Ultradent makes tips **designed to deliver** each chemistry we create. Whether you're delivering a flowable composite, a viscous gel, or thick impression material we make the perfect tip for the job. And since our tips are engineered on-site, we test each design to ensure it works perfectly with the chemistry it's intended for.



**FLEX TIP TECHNOLOGY**  
Treated for Flexibility

Our NaviTip™ tips are treated to create flexible tips with a smooth, round end, ensuring they will easily snake through canals to deliver irrigants precisely where they are needed. This treatment can create a golden hue you may see on some of our NaviTip tips.

### Check out our tips with **LOK-TITE™**

Luer Lock tips with Lok-Tite feature double threads that lock the tip into place for increased security, and wings for easy attachment and removal.



#### Black Micro™ FX™ Tip

- Accommodates various viscosities
- Flocked tip fans out to spread materials in a thin, uniform layer
- 22 ga tip

Designed for: PrimaDry, and PermaSeal.

1357 - 100pk

1434 - 500pk



#### VPS Mixing Tip

- Enables direct delivery of impression materials
- Automixing, disposable, and color coded

Designed for: Thermo Clone VPS and Chromaclone PVS.

2902 - Yellow 50pk - 2903 - Pink 50pk - 2904 - Green 50pk



#### Black Mini™ Tip

- Dispenses large volumes
- Opaque plastic preserves the flow of light-cured materials

Designed for: Ultra-Blend plus, Ultradent LC Block-Out Resin, PermaFlo, TriAway Adapter, UltraTemp, Opalescence Boost, Ultradent Diamond Polish Mint, OpalDam, Opalescence Endo, and OraSeal Caulking.

514 - 100pk

1433 - 500pk



#### Intraoral Impression Tip

- Allows precise placement
- Attaches to VPS Mixing tip

Designed for: Thermo Clone VPS and Chromaclone PVS.

2906 - 50pk




#### Black Mini™ Brush Tip

- Precise, controlled delivery of aqueous materials
- Tight, adjustable brush fibers minimize bubbles
- Unique to Ultradent

Designed for: Consepsis, Peak SE Primer, Peak-ZM, Seek/Sable Seek, Ultradent Silane, and Ultradent Universal Dentin Sealant.

1169 - 100pk

1432 - 500pk

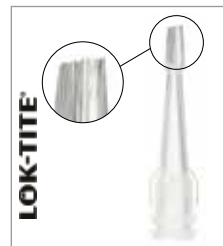

#### SoftEZ™ Tip

- Tip fibers provide visible, controlled delivery
- Brush fibers facilitate smooth application

Designed for: Enamelast.

4712 - 50pk

4711 - 100pk



#### Black Micro™ Tip Blue Micro™ Tip

- Pinpoint precision
- Narrow cannulae accurately deliver materials

Black Micro designed for: Ultra-Blend plus and EtchArrest.  
Blue Micro designed for: Ultra-Etch, Ultracare, and PrimaDry.

1085 - Black 100pk  
127 - Blue 100pk

1435 - Black 500pk  
1436 - Blue 500pk



#### White Mac™ Tip White Mini™ Tip

- Dispense large volumes
- All-plastic delivery tips

Designed for: OraSeal Caulking, OraSeal Putty, Consepsis Scrub, and Opalustre.

1361 - 100pk  
White Mac

1480 - 500pk  
White Mac

1247 - 20pk  
White Mini Laser



#### Inspiral™ Brush Tip

- Delivers viscous or filled materials smoothly via an internal helical channel and ridge
- Tight, adjustable brush fibers minimize bubbles

Designed for: Composite Wetting Resin, Peak Universal Bond, Ultradent Porcelain Etch, Ultra-Etch, UltraSeal XT plus, and UltraSeal XT hydro.

123 - 100pk

1033-I - 500pk



#### Micro 20 ga Tip

- Large-gauge cannula enables consistent flow
- Standard flowable composite delivery tip

Designed for: MTAFlow, PermaFlo, PermaFlo Purple, OpalDam, and OpalDam Green.

1252 - 100pk

1437 - 500pk



#### Metal Dento-Infusor™ Tip

- Places hemostatic agents precisely, and effectively removes superficial coagulum
- Blunt, bent cannula with padded brush enables gentle pressure on the sulcus
- Ultradent's first tip, the "MDI" remains paramount for successful tissue management

Designed for: Astringedent, Astringedent X, ViscoStat, ViscoStat Clear, and Peak Universal Bond.

2559 - 100pk

2560 - 500pk



#### Intraoral Tip

- Allows precise placement
- Attaches to dual-barrel mixing tips

Designed for: PermaFlo.

5922 - 20pk



#### Micro Capillary™ Tips

- Bright colors are easily identified against soft tissues
- The world's smallest molded tip (inside diameter: 0.008")

Designed for: Ultracare, perio materials, endodontics, and Luer Vacuum Adapter.

1120 - 5 mm 20pk

1121 - 10 mm 20pk

1802 - 5 mm 100pk



#### Micro 20 ga FX™ Tip

- Standard flowable delivery tip
- Flocked fibers fan out to evenly spread material over the entire surface

Designed for: Opalescence Boost.

1405 - 100pk

1406 - 500pk



#### SST™ (Surgical Suction Tip)

- Ideal for delicate surgeries
- Large-diameter tip opening

Designed for: Luer Vacuum Adapter for small periodontic or endodontic surgeries and controlled suction of Opalescence Boost.

1248 - 20pk



#### Blue Mini™ Dento-Infusor™ Tip

- Offers the same tissue management benefits as the Metal Dento-Infusor tip
- Allows controlled flow of drop-sized quantities
- 1.20 mm tip

Designed for: Astringedent, Astringedent X, and Consepis.

1086 - 100pk

1440 - 500pk



#### Ultradent™ Mixing Tip

- Mixes and delivers in one action

Designed for: UltraTemp REZ, UltraTemp, PermaFlo DC, and EndoREZ.

5920 - 20pk



#### Capillary Tips

- Navigate canals and substantially minimize use of paper points
- Narrow, flexible taper accesses curved canals
- 0.014" opening, 0.36 mm tip
- 0.019" opening, 0.48 mm tip

Designed for: Luer Vacuum Adapter for moisture removal from endodontic canals.

341 - 0.014" 20pk  
186 - 0.019" 20pk

3099 - 0.014" 50pk  
1425 - 0.019" 50pk